


Gamme	High-Power Diode Modules
Title	MaxxPower 405 nm
Reference	LBX-405
Description	The MaxxPower LBX-405 is a compact 405 nm laser source that exhibits the highest optical power over an industry-standard 100×40 mm footprint. Its robust design, complete with readily available safety features makes it well-suited for both laboratory and industrial use. This model is a trusted solution for polymer-curing, lithography or wide-field microscopy. Key features • ±1.0% long-term power stability; low optical noise • Modulation up to 1 MHz • Integrated control electronics Fully compatible with Oxxius MixxWave combiners for multi-mode fiber delivery in conjunction with other laser lines.

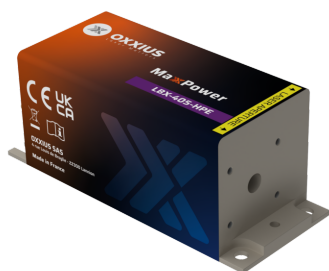

Product Variations

Part Number	Puissance
LBX-405-900-HPE-PPA	900mW
LBX-405-1200-HPE-PPA	1200mW

Optical Characteristics

Emission wavelength	405 nm
Tolerance	(±5nm tolerance)
Control modes	Automated power control (APC), Automated current control (ACC)
Optical noise	? 0.2% rms, 10Hz-20MHz bandwidth
Spectral linewidth	? 2nm
Analog modulation function - Rise / fall time	? 200ns
Analog modulation function - Bandwidth	? 1MHz
Beam diameter, 1/e ² level, 50mm from the beam aperture	2.6 x 1.2mm
Beam divergence, 1/e ² level full beam, far field	0.2 x 2.6mrad

General Specifications

Supply voltage	9V to 12V
Operating temperature	10°C to 50°C
Storage temperature and humidity	0°C to 60°C
Warm up time	?2 minutes
Interfaces	USB, RS-232, direct modulation inputs
Dimensions (laser head)	100x40x40 mm
Weight (laser head)	?330 g
Controler dimensions	106x126x56 mm
Compliance	CE (incl. IEC 60825-1) and FDA 21 CFR 1040.10 / 1040.11
Laser class	4
Warranty	24 months or 5000 hours, whichever occurs first

Options

Option 1	Multimode fiber coupling
Option 2	Heat management
Option 3	Band-pass filter
Option 4	OEM version