


FlexxRay

By Oxxius

Product range	CW & Modulated Lasers
Product name	FlexxRay 830 nm
Reference	LBX-830
Description	The FlexxRay LBX-830 is a compact 830nm laser source delivering one of the highest power available on the market over a beam of superior quality, and featuring a high-speed modulation. Its integrated control electronics and robust design suit both laboratory and industrial use. Performance, stability, and versatility make it a trusted solution for fluorescence-based microscopy, flow cytometry and illumination. Key features • TEM₀₀ spatial mode • ±0.5% long-term power stability; low optical noise • Modulation up to 150 MHz • Integrated control electronics • 100×40 mm industry-standard footprint Fully compatible with Oxxius MixxWave combiners for single-mode fiber delivery alongside other laser lines.

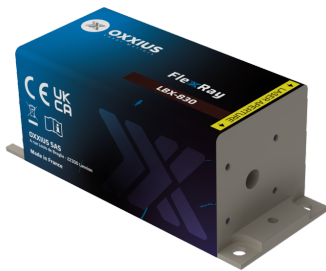

Product Variations

Part Number	Power
LBX-830-100-CSB-PPA	100mW

Optical Characteristics

Emission wavelength	830 nm
Tolerance	(±5nm tolerance)
Control modes	Automated power control, Automated current control
Optical noise	≤ 0.2% rms, 10Hz-20MHz bandwidth
Spectral linewidth	≤ 1.5nm
Digital modulation function - Rise / fall time	≤ 5ns
Digital modulation function - Bandwidth	150MHz
Analog modulation function - Rise / fall time	≤ 200ns
Analog modulation function - Bandwidth	≥ 2.5MHz
Beam diameter, 1/e² level, 50mm from the beam aperture	0.9 ±0.2mm
Beam divergence, 1/e² level full beam, far field	≤ 2.2mrad
Beam quality factor M²	≤ 1.25
Beam circularity, far field	≥ 90%
Polarization state	Linear, vertical, extinction ratio ≥ 17dB

General Specifications

Power consumption	1 W typ. / 10 W max.
Supply voltage	5V to 12V
Operating temperature	10°C to 50°C
Storage temperature and humidity	0°C to 60°C
Warm up time	≤2 minutes
Interfaces	USB, RS-232, direct modulation inputs
Dimensions (laser head)	100x40x40 mm
Weight (laser head)	≤330 g
Controller dimensions	106x126x56 mm
Compliance	CE (incl. IEC 60825-1) and FDA 21 CFR 1040.10 / 1040.11
Laser class	3B
Warranty	24 months or 5000 hours, whichever occurs first

Options

Single-mode fiber coupling
Multimode fiber coupling
Heat management
Electromechanical shutter
OEM version

Oxxius has a policy of continuous product improvement. Specifications may change without notice.