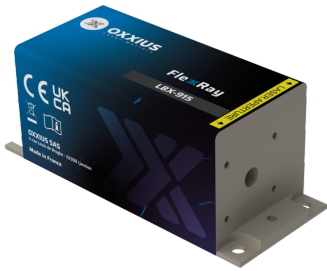


FlexxRay

By Oxxius

Product range	CW & Modulated Lasers
Product name	FlexxRay 915 nm
Reference	LBX-915
Description	The FlexxRay LBX-915 is a compact 915 nm laser source delivering one of the highest power available on the market over a beam of superior quality, and featuring a high-speed modulation. Its integrated control electronics and robust design suit both laboratory and industrial use. Performance, stability, and versatility make it a trusted solution for fluorescence-based microscopy, flow cytometry and photoluminescence.  Key features • TEM₀₀ spatial mode • ± 1 % long-term power stability; low optical noise • Modulation up to 150 MHz • Integrated control electronics • 100×40 mm industry-standard footprint Fully compatible with Oxxius MixxWave combiners for single-mode fiber delivery alongside other laser lines.

Product Variations

Part Number	Power
LBX-915-200-CSB-PPA	200 mW

Optical Characteristics

Emission wavelength	915 nm
Tolerance	(± 10 nm tolerance)
Power stability over 8 hours, temperature within $\pm 3^{\circ}\text{C}$	± 1 %
Control modes	Automated power control, Automated current control
Optical noise	≤ 0.2 % rms, 10 Hz-20 MHz bandwidth
Spectral linewidth	≤ 1.2 nm
Digital modulation function - Rise / fall time	≤ 5 ns
Digital modulation function - Bandwidth	150 MHz
Analog modulation function - Rise / fall time	≤ 200 ns
Analog modulation function - Bandwidth	≥ 2 MHz
Beam diameter, $1/e^2$ level, 50mm from the beam aperture	0.8 ± 0.2 mm
Beam divergence, $1/e^2$ level full beam, far field	≤ 2 mrad
Beam quality factor M^2	≤ 1.25
Beam circularity, far field	≥ 90 %
Polarization state	Linear, vertical, extinction ratio ≥ 17 dB

General Specifications

Power consumption	10 W max.
Supply voltage	5 V to 12 V
Operating temperature	10°C to 50°C
Storage temperature and humidity	0°C to 60°C
Warm up time	≤ 2 minutes
Interfaces	USB, RS-232, direct modulation inputs
Dimensions (laser head)	100x40x40 mm
Weight (laser head)	≤ 330 g
Controller dimensions	106x126x56 mm
Compliance	CE (incl. IEC 60825-1) and FDA 21 CFR 1040.10 / 1040.11
Laser class	3B
Warranty	24 months or 5000 hours, whichever occurs first

Options

Single-mode fiber coupling
Multimode fiber coupling
Heat management
Electromechanical shutter
OEM version

Oxxius has a policy of continuous product improvement. Specifications may change without notice.