


MaxxPower

By Oxxius

Product range	High-Power Diode Modules
Product name	MaxxPower 980 nm
Reference	LBX-980
Description	The MaxxPower LBX-980 is a compact 980 nm laser source that exhibits the highest optical power over an industry-standard 100×40 mm footprint. Its robust design, complete with readily available safety features makes it well-suited for both laboratory and industrial use. This model is a trusted solution for illumination or optical pumping. Key features • $\pm 1.0\%$ long-term power stability; low optical noise • Modulation up to 1 MHz • Integrated control electronics Fully compatible with Oxxius MixxWave combiners for multi-mode fiber delivery in conjunction with other laser lines.

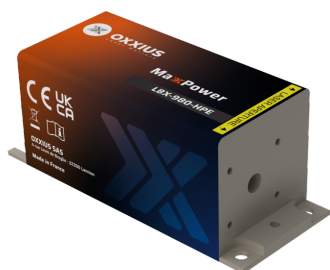

Product Variations

Part Number	Power
LBX-980-1000-HPE-PPA	1000mW

Optical Characteristics

Emission wavelength	980 nm
Tolerance	(±5nm tolerance)
Control modes	Automated power control (APC), Automated current control (ACC)
Optical noise	≤ 0.8% rms, 10Hz-20MHz bandwidth
Spectral linewidth	≤ 2nm
Analog modulation function - Rise / fall time	≤ 200ns
Analog modulation function - Bandwidth	≥ 1MHz
Beam diameter, 1/e² level, 50mm from the beam aperture	3.5 x 1.5mm
Beam divergence, 1/e² level full beam, far field	0.3 x 12mrad

General Specifications

Supply voltage	5V to 12V
Operating temperature	10°C to 50°C
Storage temperature and humidity	0°C to 60°C
Warm up time	≤2 minutes
Interfaces	USB, RS-232, direct modulation inputs
Dimensions (laser head)	100x40x40 mm
Weight (laser head)	≤330 g
Controler dimensions	106x126x56 mm
Compliance	CE (incl. IEC 60825-1) and FDA 21 CFR 1040.10 / 1040.11
Laser class	4
Warranty	24 months or 5000 hours, whichever occurs first

Options

Multimode fiber coupling
Heat management
Band-pass filter
OEM version

Oxxius has a policy of continuous product improvement. Specifications may change without notice.